

Installing the Onscreen Keyboard App On Your T10 Device

The Onscreen Keyboard application is a stand-alone Android application from DynaVox that gives the T10 user the ability to replace the standard Android onscreen keyboard with custom keyboards based on DynaVox page content. This permits the use of familiar access methods, word prediction, and user interface layouts to insert text in non-DynaVox Android applications.

The Onscreen Keyboard version 1.1 application is now available at myDynaVox.com.

**NOTE**


You must have Compass 1.3 installed on your T10 in order to use Onscreen Keyboard 1.1

You can download the Onscreen Keyboard version 1.1 app directly to your T10 device, then install it by following the instructions in Section I of this document.


Alternatively, you can download the Onscreen Keyboard 1.1 app to your PC, then install to your T10 device following the instructions in Section II. (*II. Downloading to Your PC and Installing on Your T10*).

I. Downloading to the T10


1. On your T-10 device, go to myDynaVox.com.
2. Log in using your myDynaVox username and password.
3. Select **Onscreen Keyboard** to start downloading the Onscreen Keyboard app.

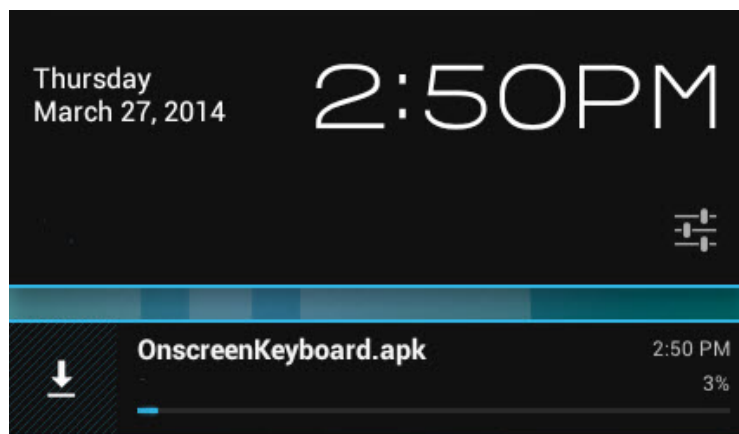


4. Select **Browser**.

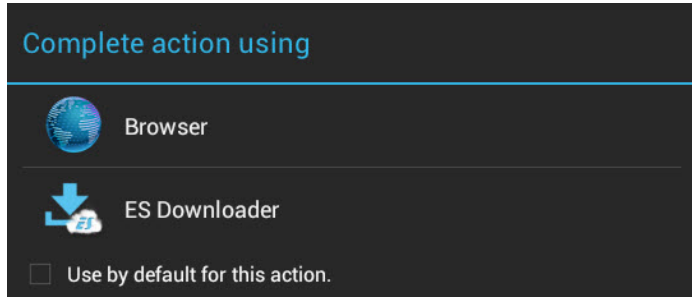


5. Wait for the download to complete.

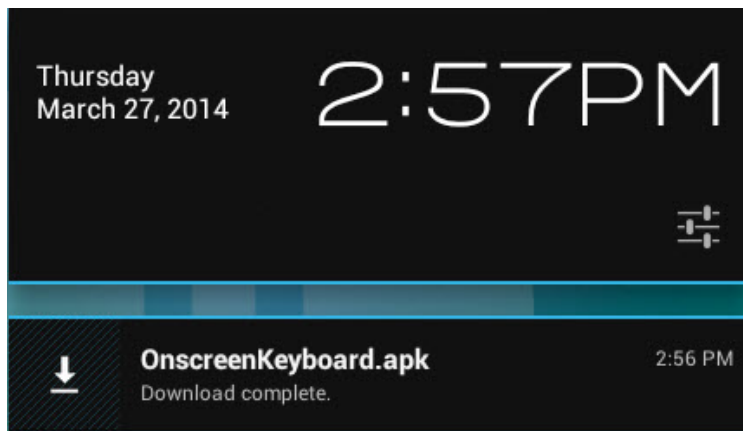


6. Open the *Settings* menu on your T-10 desktop. 

7. Select **Storage** on the left..

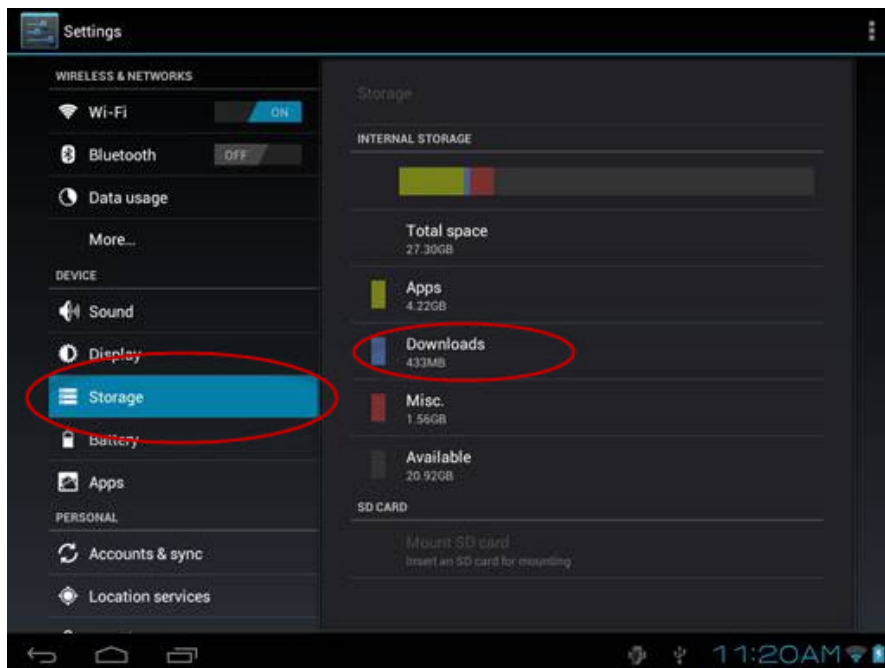


8. Select **Downloads** on the right.
9. Select **Onscreen Keyboard.apk**.
10. On the *Replace application* dialog that opens, select **OK**.

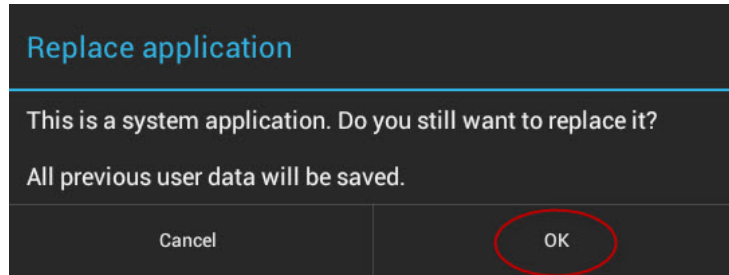

The *Onscreen Keyboard* screen will open.

11. Select **Install** on the *Onscreen Keyboard* screen

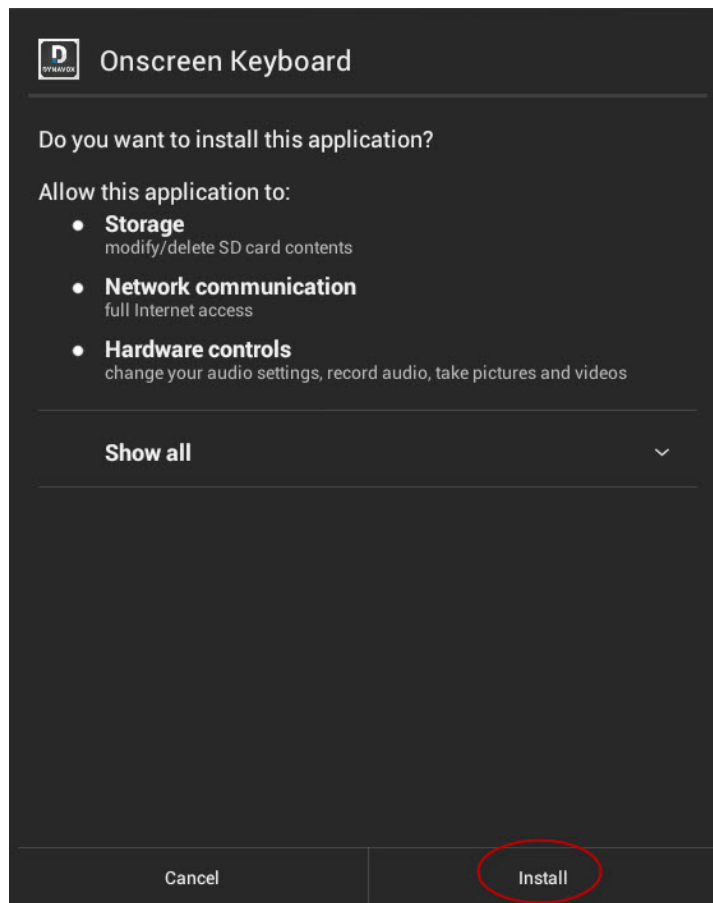

12. Select **Done** on the *Complete the Installation* dialog.

Onscreen Keyboard 1.1 is now installed on your T-10 device.

II. Downloading to Your PC and Installing on Your T10


1. On your PC web browser, navigate to *myDynaVox.com*.
2. Log in using your myDynaVox username and password.
3. Select **Onscreen Keyboard 1.1** to start downloading the Onscreen Keyboard app.




NOTE


IMPORTANT! Save the file either to your Desktop or to a folder where you can easily find it.


4. Connect your T10 to your PC, using the USB to micro USB cable that was shipped with your device.
5. On your PC, using Windows Explorer, copy the *OnscreenKeyboard.apk* file onto your T10:
 - a. Select the T10 icon in the left column, under *Computer*.
 - b. Open the *Internal Storage* folder.
 - c. Open the *Download* folder.

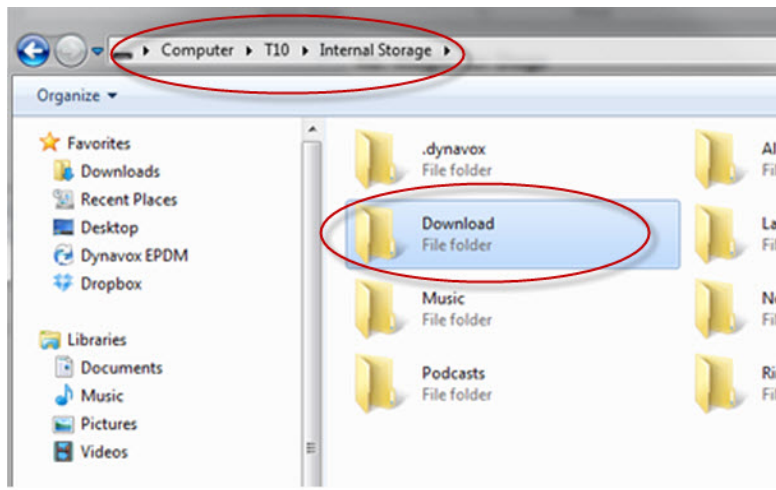


- d. Copy the *.apk* file into the *Download* folder.

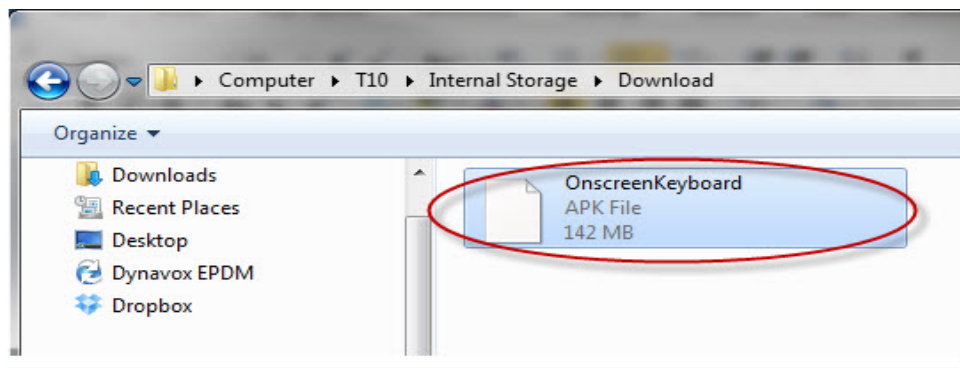


6. On your T10, go to the app browser. ( .)
7. Open *ES File Explorer* .
8. Open the *Download* folder.
9. Select **Onscreen Keyboard.apk**.

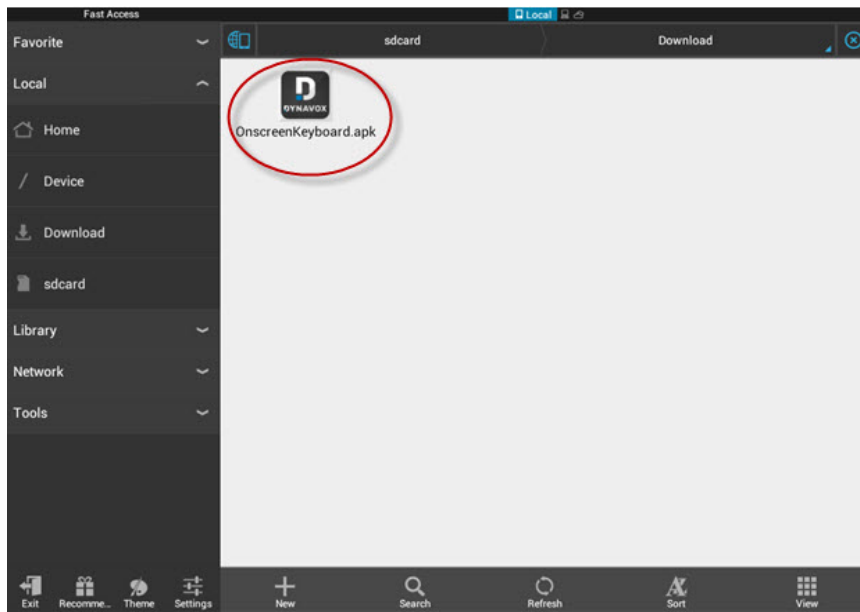

10. Select **OK** on the *Replace application* dialog that opens.

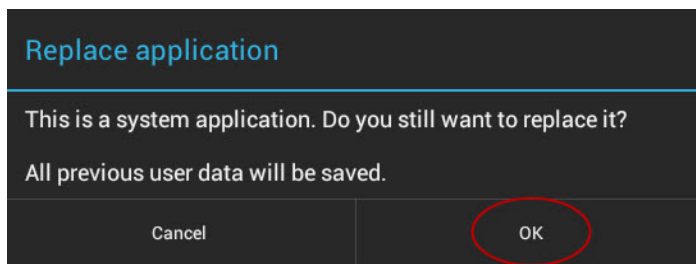

11. Select **Install** on the *Onscreen Keyboard* dialog that opens.

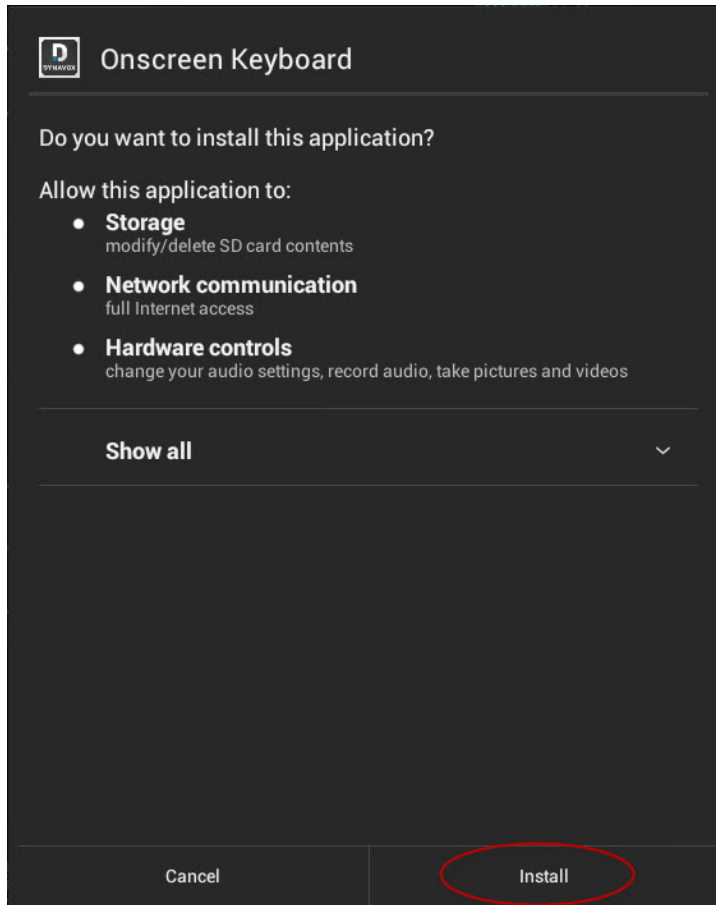

12. Select **Done** on the *Complete the Installation* dialog.

Onscreen Keyboard 1.1 is now installed.