


AccessIT 3

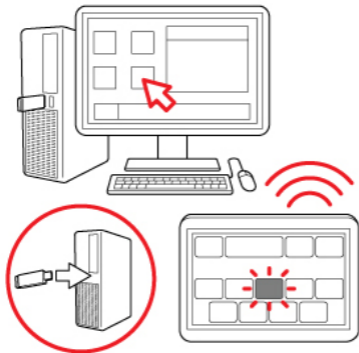

This document, as well as installation instructions, available online in the following languages: da-DK, de-DE, en-US, es-MX, fr-FR, nb-NO, nl-NL, sv-SE, zh-CN. Visit

<http://www.tobii dynamox.com/AccessIT>

Safety

The AccessIT 3 has been tested and approved as compliant to all the Specifications and Standards listed in the

Copyright ©Tobii AB (publ). Not all products and services offered in each local market. Specifications are subject to change without prior notice. All trademarks are the property of their respective owners.

Compliance Information section. To ensure safe operation, there are a few safety warnings to bear in mind:



No modification of this equipment is allowed.

Do not open the AccessIT 3. Non-compliance will result in loss of Warranty! There are no user-serviceable components inside.



The AccessIT 3 device shall not be exposed to or used in rain or wet conditions.



The AccessIT 3 is not to be used as a life supporting Device, and it shall not be relied on in case of loss of function due to power loss or other causes.



Do not rely on the device for emergency calls or banking transactions. We recommend having multiple ways to communicate in emergency situations. Banking transactions should only be carried out with a system recommended by, and approved according to the standards of, your bank.



Young children or people with cognitive disabilities should not have access to, or the use of, the AccessIT 3 device without parental or guardian supervision.



The AccessIT 3 is a compact electronic device. As such, it is composed of numerous separate, assembled parts. In the hands of a child, the AccessIT 3 or its parts could constitute a choking hazard or other danger to the child.



Any use of the AccessIT 3 outside the intended use and together with any third party software or hardware that changes the intended use is a risk and nothing that Tobii Dynavox can take any responsibility for.

Compliance Information



The AccessIT 3 is CE-marked, indicating compliance with the essential health and safety requirements set out in European Directives.

FCC Statement

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.



Modifications not expressly approved by Tobii Dynavox could void the user's authority to operate the equipment under FCC rules.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be

determined by disconnecting and reconnecting the AccessIT 3, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into a wall outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or experienced radio/TV technician for assistance.

Industry Canada Statement

This Class B digital apparatus complies with Canadian ICES-003.

Cet appareil numérique de la classe B est conforme à la norme NMB-003 du Canada.

CE Statement & Standards

The AccessIT 3 complies with the following Directives:

- 2011/65/EU (RoHS 2) - Restriction of Hazardous Substances Directive
- 2006/95/EC (LVD) - Low Voltage Directive
- Radiocommunication Standard 2008 (ACMA), AS/NZS CISPR 32:2013

The AccessIT 3 complies with the following standards:





- FCC, CFR Title 47, Part 15 Subpart B: 2015, CLASS B
- ICES-003 Issue 6:2016 Class B
- EN 55032: 2012+AC: 2013 Class B
- EN 55024:2010+A1:2015
- EN 61000-3-2: 2014
- EN 61000-3-3: 2013

Disposing of the AccessIT 3

Do not dispose of the AccessIT 3 in general household or office waste. Follow your local regulations for the disposal of electrical and electronic equipment.

Symbols and Markings

This section provides information about the symbols that are used on the AccessIT 3 or its packaging.

Symbol or Marking	Description
	Dispose of in accordance with your country's requirements.
	CE is the abbreviation of the European Communities and this mark tells customs officials in the European Union that the product complies with one or more of the EC Directives.
	This mark is a certification mark employed on electronic products manufactured or sold in the United States which certifies that the electromagnetic interference from the device is under limits approved by the Federal Communications Commission.
	Australian Electrical Equipment Safety System Regulatory Compliance Mark.

Technical Specifications

Technical Specifications	AccessIT 3
Unit Size (Length × Height × Depth)	70mm × 23mm × 9mm 2.75" × 0.9" × 0.35"
Unit Weight	14 g 0.5 oz
Active range of Bluetooth transmission	10 m 32 feet
Works with (Operating System)	Windows 10 OS X, macOS, iOS
Connectivity	USB 2.0 or higher

Customer Support

For support, please contact your local representative or Tobii Dynavox support. Find our Support pages online at: www.tobiidynavox.com or www.myTobiiDynavox.com.

Contact Information:

Tobii Dynavox
Karlsrovägen 2D
182 53 Danderyd
Sweden

+46 8 663 69 90

Tobii Dynavox LLC
2100 Wharton Street, Suite 400
Pittsburgh, PA 15203
USA

+1-800-344-1778